

# ARCTIKO

The art of simplicity®

## ULUF 750

with DUAL SINGLE COMPRESSOR TECHNOLOGY



### LIMITED OFFER!\*

This unit also uses the ARCTIKO single compressor technology, BUT the ULUF 750 has two refrigeration systems!

With the dual systems working together the pull-down time from room temperature to -80°C is 3.5 hours, versus 8 hour pull-down time for the conventional "Cascade system" you have.

All while using the same amount of electrical consumption (technically you would be saving in energy cost due to the quick recovery time of the ULUF 750).

We did a 2 minute door opening test on the ULUF 750 and it recovered back to -80°C in 35 minutes. Other freezers will take well over 2 hours!


MODEL	ULUF 750
Item Code	DAI 1416
Dimensions exterior WxDxH (mm)	1030x885x2089
Capacity (L)	680
Temperature range (°C)	-86 / -40
Power supply (V)	230
Frequency (Hz)	50 / 60

-86°C

€ 7.950,-\*

#### IN CASE OF TROUBLE

#### NORMAL OPERATION

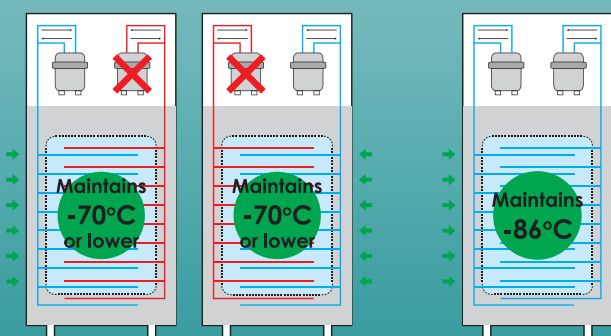

#### DOUBLE COOLING SYSTEM / Normal Operation:

Two independent Single compressor cooling systems work simultaneously. Each system is able to reach ultra-low temperature.


#### IN CASE OF TROUBLE / One System down:

One of the Systems can maintain -70°C or lower.

For more information, visit [www.arctiko.com](http://www.arctiko.com)

**bartelt**

IHR LABOR-KOMPLETTAUSSTATTER  
VERKAUF - GERÄTESERVICE - SOFTWARE

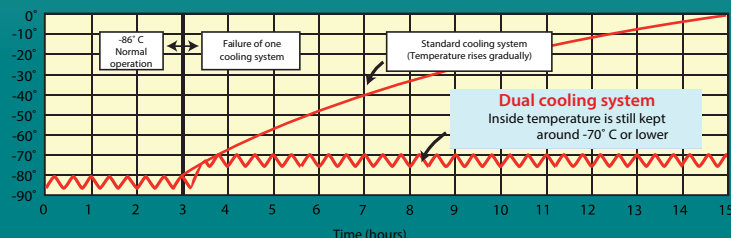

\* All Prices are Net Prices.

We reserve the right to change specifications without notice.  
Subject to confirmation, availability and errors.  
Valid until 30.09.2017 including transport.